




# Integrating Four Core Technologies TO SHAPE THE FUTURE OF MEDICINE

ASAHI® IS THE INDUSTRY LEADER IN GUIDE WIRE TECHNOLOGY. INTEGRATION OF 4 CORE PROPRIETARY TECHNOLOGIES PROVIDES PURPOSEFUL OPTIMIZATION OF GUIDE WIRE PERFORMANCE FOR BOTH ROUTINE AND COMPLEX PROCEDURES.

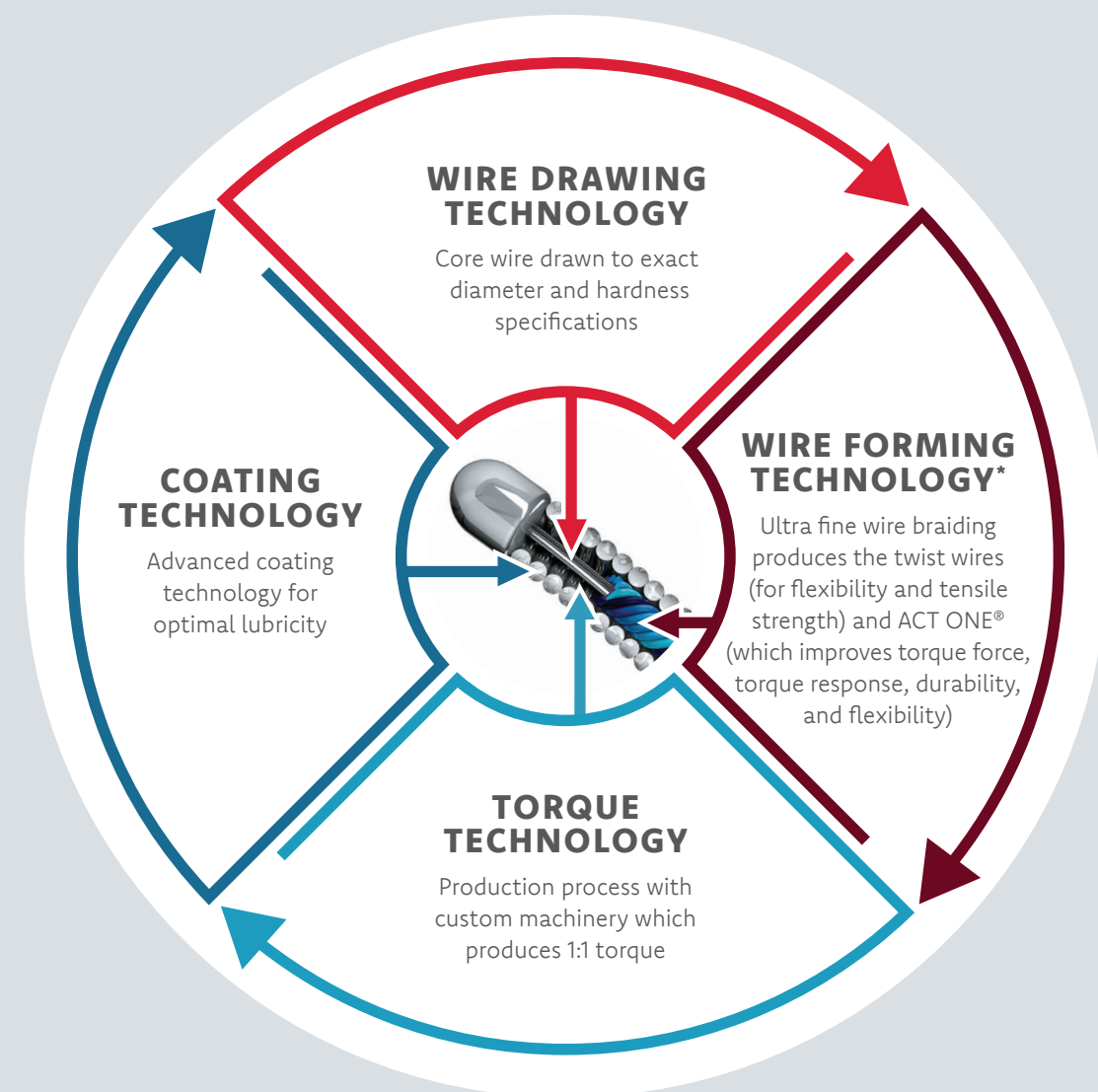

\*Available in select ASAHI guide wires

## ORDERING INFORMATION

FRONTLINE GUIDE WIRES					
PRODUCT	CATALOG NO.	TIP SHAPE	LENGTH (CM)	SPRING COIL LENGTH (CM)	COATING
ASAHI SION blue 190 cm	AHW14R104S / AHW14R104J	Straight/J	190	20	Hybrid Coating: Hydrophilic* with Hydrophobic Tip
ASAHI SION blue 300 cm	AHW14R304S / AHW14R304J	Straight/J	300		
MINAMO 190cm	AH14R026S	Straight	190	20	Hybrid Coating: Hydrophilic* with Hydrophobic Tip
	AH14R026J / AH14R026P	J/Pre-Shape			
MINAMO 300cm	AH14R326S	Straight	300	28	Hydrophilic* coating over spring coil
	AH14R326J / AH14R326P	J/Pre-Shape			
ASAHI SION	AHW14R101S	Straight	190	12	Polymer Jacket & Hydrophilic*
ASAHI SION 300 cm	AHW14R301S	Straight	300		
ASAHI SION black 190cm	APW14R010S	Straight	190	4	Hydrophobic (Silicone)
	APW14R010J / APW14R010P	J/Pre-Shape			
ASAHI SION black 300cm	APW14R310S	Straight	300		
	APW14R310J / APW14R310P	J/Pre-Shape			
GRAND SLAM	AG141002 / AG141002J	Straight/J	180		
GRAND SLAM 300 cm	AG141302 / AG141302J	Straight/J	300		

CHRONIC TOTAL OCCLUSION GUIDE WIRES					
PRODUCT	CATALOG NO.	TIP SHAPE	LENGTH (CM)	SPRING COIL LENGTH (CM)	COATING
FIELDER XT-A	APW14R009S	Straight	190	16	Polymer Jacket & Hydrophilic*
FIELDER XT-A 300 cm	APW14R309S	Straight	300		
FIELDER XT-R	APW14R005S	Straight	190	16	Polymer Jacket & Hydrophilic* Hydrophilic coating over polymer jacket: 17cm
FIELDER XT-R 300 cm	APW14R305S	Straight	300		
ULTIMATEbros 3	AHW14S003S	Straight	180	11	Hydrophilic* with Uncoated Distal Tip
ULTIMATEbros 3 300 cm	AHW14S303S	Straight	300		
CONFIANZA Pro	AGH143090	Straight	180	20	Hydrophilic* with Uncoated Distal Tip
CONFIANZA Pro 300 cm	AGH143390	Straight	300		
CONFIANZA Pro 12	AGH143091	Straight	180	20	Hydrophilic* with Uncoated Distal Tip
CONFIANZA Pro 12 300 cm	AGH143391	Straight	300		
ASAHI Gaia Next1	AH14R019P	Pre-Shape	190	15	Hydrophilic* with Uncoated Distal Tip (SLIP-COAT®) 40 cm
ASAHI Gaia Next2	AH14R020P	Pre-Shape	190		
ASAHI Gaia Next3	AH14R021P	Pre-Shape	190		

PURPOSE BUILT GUIDE WIRES					
PRODUCT	CATALOG NO.	TIP SHAPE	LENGTH (CM)	SPRING COIL LENGTH (CM)	COATING
ASAHI Gladius Mongo, 190 cm	AP14R023S / AP14R023P	Straight/Pre	190	8.5	Polymer Jacket & Hydrophilic*
ASAHI Gladius Mongo, 300 cm	AP14R323S / AP14R323P	Straight/Pre	300		
ASAHI SUOH 03	AHW14R013S / AHW14R013P	Straight/Pre	190	19	Hydrophilic*
ASAHI SUOH 03 300 cm	AHW14R313S / AHW14R313P	Straight/Pre	300		
ASAHI SION	AHW14R101S	Straight	190	28	Hydrophilic*
ASAHI SION 300 cm	AHW14R301S	Straight	300		
ASAHI RG3	AHW10S302S	Straight	330	8	Hydrophilic*
ASAHI EXTENSION 165 cm	AG149001	Straight	165	-	-



Toll Free: 855-286-9473 Fax: 949-756-8146 WWW.ASAHI-INTECC-US.COM

**Indications for Use** ASAHI PTCA Guide Wires are intended to facilitate the placement of balloon dilation catheters during percutaneous transluminal coronary angioplasty (PTCA) and percutaneous transluminal angioplasty (PTA). The ASAHI PTCA Guide Wires are not to be used in the neurovasculature.

© 2020 ASAHI INTECC CO., LTD. "ASAHI", "ACT ONE", "SION", "ASAHI Gaia", and "FIELDER" are trademarks or registered trademarks of ASAHI INTECC CO., LTD. in Japan and other countries.  
\*Coated with SLIP-COAT® coating  
Some materials / graphics were provided by Vascular Perspectives – www.vascularperspectives.com.

K21200\_B\_LC\_2.2  
MKT0394 Rev. B

# CORONARY GUIDE WIRES

MASTER THE EVERYDAY



ASAHI INTECC USA, INC.

## FRONTLINE GUIDE WIRES



**ASAHI SION® blue**

Frontline guide wire with proprietary ACT ONE construction for greater tip flexibility and support.  
(Easily-discernible shaft color for differentiation from other ASAHI wires.)

TIP LOAD	TIP RADIOPACITY	COATING
0.5 gf	3 cm	Hybrid coating Hydrophobic coating: tip to 1.5 cm Hydrophilic coating 18.5 cm

**MINAMO®**

A workhorse frontline guide wire with ACT ONE and Nitinol tip. Superior torque response with the resiliency of nitinol. Less whip than other similar guide wires.

TIP LOAD	TIP RADIOPACITY	COATING
0.5 gf	3 cm	Hybrid coating Hydrophobic coating: tip to 2.5 cm Hydrophilic coating: 50 cm

**ASAHI SION**

Frontline guide wire with ACT ONE for excellent tip durability and torque performance in tight anatomy. Flexible shaft for greater trackability.

TIP LOAD	TIP RADIOPACITY	COATING
0.7 gf	3 cm	Hydrophilic coating: 28 cm

**ASAHI SION black**

A workhorse frontline guide wire with ACT ONE for torque response and increased durability. Also has a polymer jacket for lubricity. High tip flexibility improves trackability in tortuous anatomy.


TIP LOAD	TIP RADIOPACITY	COATING
0.8 gf	3 cm	Hydrophilic coating over polymer jacket

**GRAND SLAM®**

Frontline guide wire with extra support for delivery of interventional devices.

TIP LOAD	TIP RADIOPACITY	COATING
0.7 gf	4 cm	Hydrophobic coating over spring coil

## CHRONIC TOTAL OCCLUSION GUIDE WIRES

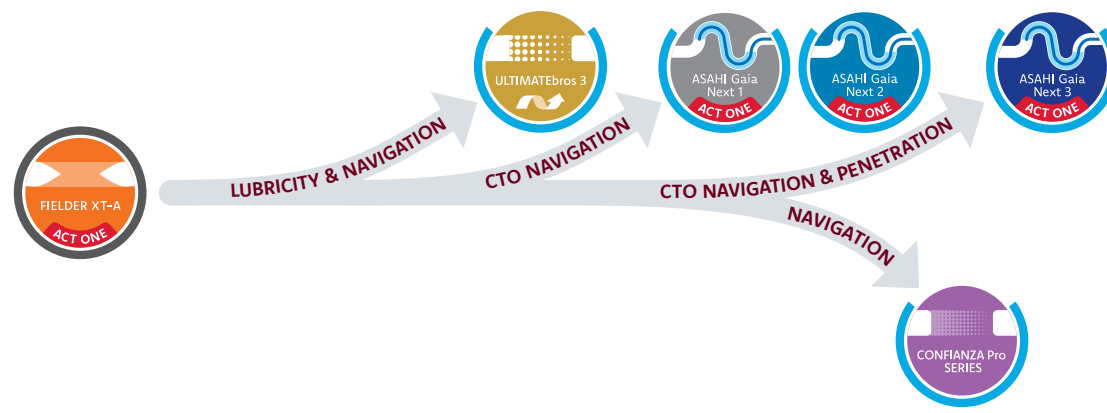

**FIELDER® XT-A** CTO

Guide wire with tapered durable tip, facilitating entry into the chronic occluded lesions.

TIP LOAD	TIP RADIOPACITY	POLYMER JACKET LENGTH	TIP OUTER DIAMETER	COATING
1.0 gf	16 cm	17 cm	0.010 inch (0.26 mm)	Hydrophilic coating over polymer jacket

**FIELDER XT-R** CTO

Tapered guide wire with a flexible tip and lubricious polymer coating.

TIP LOAD	TIP RADIOPACITY	POLYMER JACKET LENGTH	TIP OUTER DIAMETER	COATING
0.6 gf	16 cm	17 cm	0.010 inch (0.26 mm)	Hydrophilic coating over polymer jacket

**ULTIMATEBros® 3** CTO

Long hydrophilic coating maintains high maneuverability, allowing for improved wire manipulation in highly stenosed lesions. Fine shaping improves vessel selectivity and reduces the risk of false lumen expansion.

TIP LOAD	TIP RADIOPACITY	COATING
3.0 gf	11 cm	Hydrophilic 40 cm - uncoated distal tip

**CONFIANZA® Pro Series** CTO

Tapered guide wire for crossing complex lesion with heavy calcification or dense fibrous tissues.


TIP LOAD	TIP RADIOPACITY	TIP OUTER DIAMETER	COATING
9.0 / 12.0 gf	20 cm	Pro, Pro12 0.009 inch (0.23 mm)	Hydrophilic 20 cm - uncoated distal tip

**ASAHI Gaia® Next Series (One, Two, Three)** CTO

Guide wire with composite core and XTRAND coil to improve navigation and increase anti-trapping power.

TIP LOAD	TIP RADIOPACITY	TIP OUTER DIAMETER	COATING
One: 2.0 gf Two: 4.0 gf Three: 6.0 gf	15 cm	One: 0.011 inch (0.27mm) Two: 0.012 inch (0.30mm) Three: 0.012 inch (0.30mm)	Hydrophilic 40 cm - uncoated distal tip

## PURPOSE BUILT GUIDE WIRES

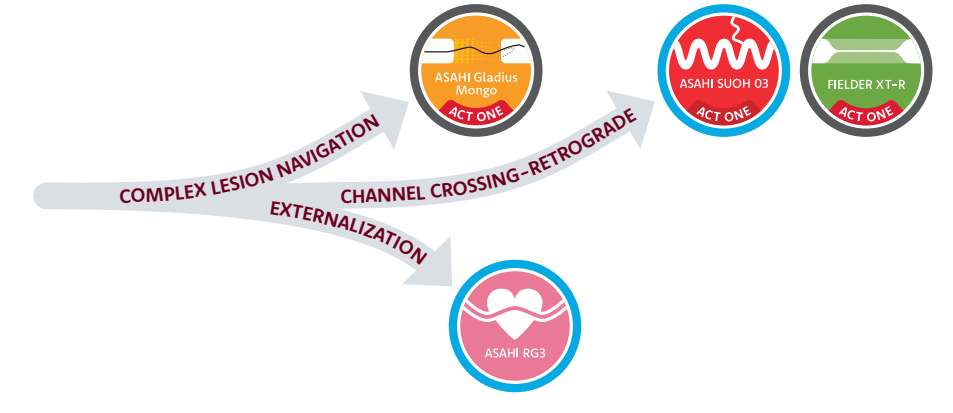

**ASAHI Gladius® Mongo®**

A new guide wire designed for maximized steerability as well as providing the added ability to proplase into a small controlled area when working in the subintimal space. ES provides Extra Support for improved push force and support.

TIP LOAD	TIP RADIOPACITY	COATING
3.0 gf	3 cm	Hydrophilic coating over polymer jacket

**ASAHI SUOH® 03**

Flexible tip and shaft provide superb trackability and crossability in severe tortuosity.

TIP LOAD	TIP RADIOPACITY	COATING
0.3 gf	3 cm	Hydrophilic coating: 52 cm

**ASAHI RG3**


For guide wire externalization only.

TIP LOAD	TIP RADIOPACITY	SHAFT DIAMETER	GUIDE WIRE LENGTH	COATING
3.0 gf	3 cm	0.010 inch (0.26 mm)	330 cm	Hydrophilic coating distal 170 cm Hydrophobic coating proximal 160 cm

**ASAHI EXTENSION® 165**

Designed to add 165 cm to the end of an ASAHI non-300cm wire.

### WIRE COATINGS LEGEND

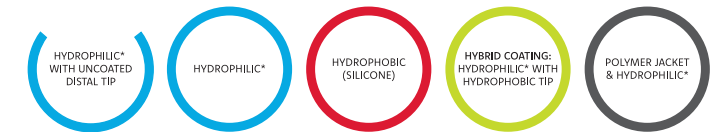

CTO Indicates guide wires with a Chronic Total Occlusion Indication